## User's Manual

### Pointer for Touch Panel Color LCD Monitor

Thank you very much for choosing an EIZO product.

Please read this User's Manual carefully to familiarize yourself with safe and effective usage procedures.

Please retain this manual for future reference.

# Package Contents

Please contact your local dealer for assistance if any of the listed items are missing or damaged.

- Touch Panel Pointer
- Pointer Holder
- Screw for pointer installation × 1
- User's Manual

# How to install

A driver is required to install the touch panel pointer and pointer holder.


### **Touch Panel Pointer installation**


Attach the wire lead of the pointer to the monitor by using the supplied screw.


- When the monitor has the stand-mounting holes on the rear, attach the TP1 into the lower hole.
- When the monitor has the hole exclusively for TP1, attach the TP1 into it. (About the location of hole for TP1, refer to User's Manual of the monitor.)

#### Pointer Holder installation

- (1) Remove the upper screw from the rear of the monitor.
- (2) Insert the stopper of the holder into the ventilation slot of the monitor.
- (3) Fix the holder on the monitor by using the screw removed at step (1).
- (4) House the pointer into the holder.

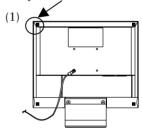

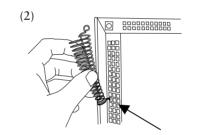

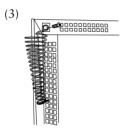



# **Specifications**

	Dimensions	Weight
Touch Panel Pointer	140 mm (L)× φ 10 mm (5.5" (L) × φ 0.39")	50 g (1.76 oz)
Pointer Holder	24.1 mm (W) × 86.3 mm (H) × 28 mm (D) (0.95" (W) × 3.40" (H) × 1.10" (D))	10 g (0.35 oz)

#### ■ Dimensions

Touch Panel Pointer

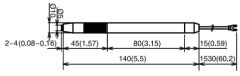

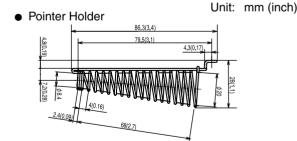

TP1 is designed for pointing on the screen. The drawing touch and dragging it on the touch panel is not guaranteed.

## **EIZO** Corporation

153 Shimokashiwano, Hakusan, Ishikawa 924-8566 Japan Phone: +81 76 277 6792 Fax: +81 76 277 6793

# EIZO Inc.

5710 Warland Drive, Cypress, CA 90630 Phone: +1 562 431 5011 Fax: +1 562 431 4811

### EIZ口株式会社

〒924-8566 石川県白山市下柏野町153番地

http://www.eizo.com http://www.eizo.co.jp



6th Edition-January, 2014

Printed in Japan.