


# Stand Mounting Instructions for the RadiForce® LX470W

The stand is designed for the RadiForce® LX470W. Please ensure that all monitor weight and stability requirements are met before using the stand with other models. (VESA mount 400mm x 200mm)

Order No.: 6GF69888BA03


Description: FST4700

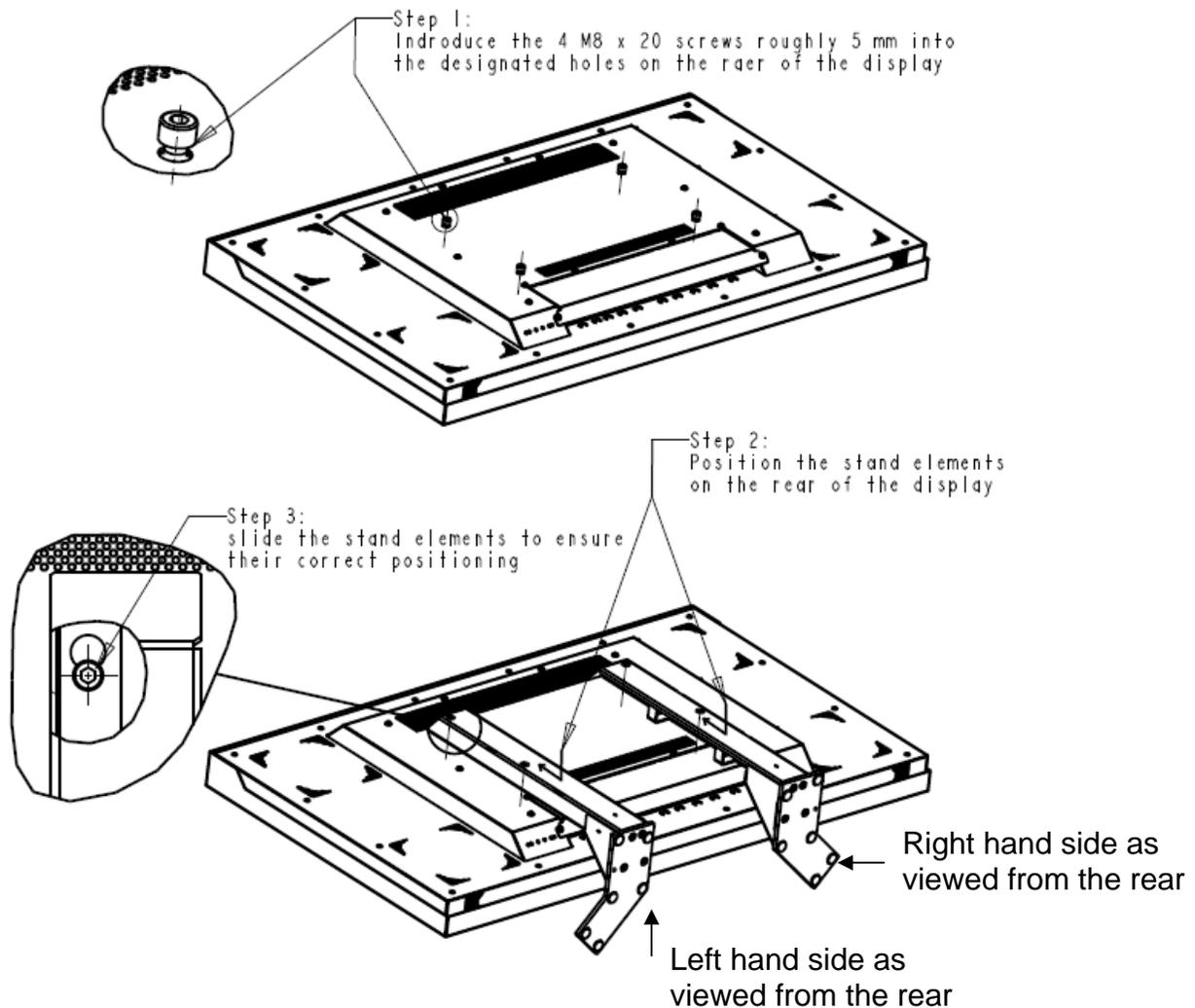

Tighten the 4 screws with an Allen wrench (18 Nm)

## CAUTION

- The device must be removed from the packaging and carried by at least 2 persons.
- Wear appropriate protection to prevent injuries should the device be dropped.

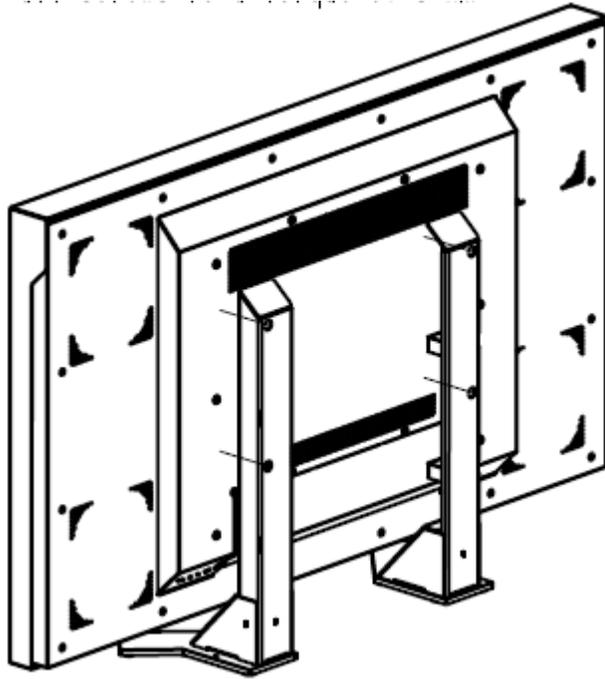

Place the monitor in its final intended location.